

Tubulicid Plus Endo[®]

– FOR CLEANING OF ROOT CANALS AND
DENTINAL TUBULE APERTURES


Tubulicid Plus Endo is specially formulated for routine and dedicated root canal therapy, combining a selective demineralizing effect with a surface active treatment.

For effective cleaning of root canals it is important that the smear layers produced during mechanical root canal preparation are removed and that the debris blocking the apertures of dentinal tubules is cleared away.

With its combination of low surface tension and EDTA (3%), Tubulicid Plus Endo facilitates penetration of the dentinal tubules which effectively removes both the outer continuous smear layers and smear plugs, increasing the adhesion of the filling material.

Facts and benefits

- **Root canal cleanser** - Will clean canals and dentinal tubules with EDTA (3%)
- **Antibacterial agent** - Will disinfect infected root canals
- **Removes smear layer** - Clears debris caused by mechanical root canal preparation
- **Low surface tension** - Facilitates penetration of dentinal tubules

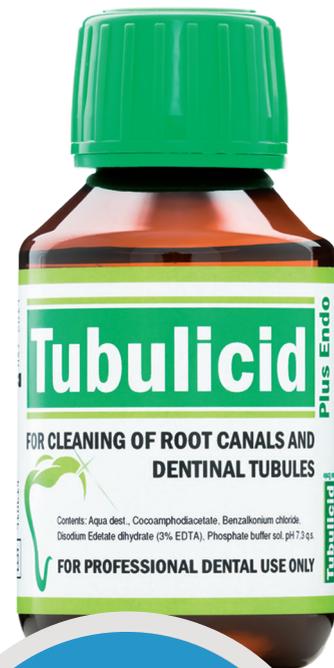

Formulated during years of research by the “father of dentinal tubules”
Dr. Martin Brännström.

 **Dental Therapeutics**

phone: +46(0)8 716 81 41

email: info@therapeutics.se

web: www.therapeutics.se